

Test_Facility_Robot_UK

A healthcare professional in blue scrubs is shown from the chest down, holding a silver smartphone in their right hand and a black pen over a spiral-bound notebook in their left hand. A stethoscope is draped around their neck. The background is blurred, showing what appears to be a window or outdoor area. Overlaid on the image are numerous white circular icons connected by thin white lines, forming a network-like diagram. These icons represent various aspects of healthcare and technology: a doctor's head, a first aid kit, a person at a computer, a bar chart, a person wearing a headset, a magnifying glass, a gear, crossed wrench and screwdriver, a document with a checkmark, a mobile phone, a factory/building, a person sitting at a desk, a calculator, a document with a checkmark, a person walking, a truck with a cross, a laptop, and a central clock face. At the bottom center, there is a stylized blue lightning bolt icon.



Executive Summary

BeaconMedaes Service Technician Robot_UK Tester tested the Medical Gas Pipeline System indicated by this report in accordance with HTM02-01 Part B section 10 (Maintenance) and/or accordance with the facilities own specification (where requested) and other applicable standards or agencies. The tests performed on this system were found to be within test specifications unless otherwise noted. The results, dates and locations of the tests performed and specific data on each test can be found on the Test Summary pages.

The testing of this system may include all or part of the tests offered by BeaconMedaes. Only the tests indicated by the summary and data pages were performed on the dates indicated. The results of the tests reported represent BeaconMedaes findings on the dates and locations indicated by the report forms only. System components include:

1. Air System : 1 pcs

Inspection Report

Test Date: 28-Jan-25

Air System: Oil Lub Screw: Medical: 7.5 Hp / 5.5 KW: Triplex: BeaconMedaes: 4233600592: mAL
HTM02-01 AC GA5 dMED35, s/n MNL_2498

Location: Building 1 : Floor 1

01. General

Data Point	Result
Plant Room Condition	Pass
Plant access	Pass
Plant fixing and bonding	Pass
Compressor manufacturer	Atlas Copco
Dryer manufacturer	BeaconMedaes
Relevant staff advised of system test and likely alarms	Pass
ERM is online and able to supply air to the MGPS	Pass
Safety valve(s) in date	Pass
Receiver(s) condition is good	Pass
Dryer condition is good and operates correctly	Pass
Dryer line pressure (bar)	0.000
Gauge(s) operate correctly	Pass
Condition of electrical components are acceptable	Pass
Isolate each compressor and check running conditions are acceptable	Pass
Compressor(s) are free from oil leaks	Pass
Alarms and indications operate correctly	Pass
Plant run hours	0.000
Pressure switch fault level (bar)	0.000
Backup pressure switch operates correctly	Pass
Filters are clean and in date	Pass
Filter auto drains are operating correctly	Pass
Equipment status is online	Pass
Atmospheric dew point level	0.000

02.1 Compressor 1

Data Point	Result
Compressor run hours	0.000
Cut in pressure (bar)	0.000
Cut out pressure (bar)	0.000
Running current loaded (A)	0.000
Running current unloaded (A)	0.000
Oil level check	Pass
Oil temperature	0.000

02.2 Compressor 2

Data Point	Result
Compressor run hours	0.000
Cut in pressure (bar)	0.000
Cut out pressure (bar)	0.000
Running current loaded (A)	0.000
Running current unloaded (A)	0.000

Oil level check	Pass
Oil temperature	0.000

02.3 Compressor 3

Data Point	Result
Compressor run hours	0.000
Cut in pressure (bar)	0.000
Cut out pressure (bar)	0.000
Running current loaded (A)	0.000
Running current unloaded (A)	0.000
Oil level check	Pass
Oil temperature	0.000



03. Dynamic Risk Assessment

Data Point	Result
Risk Assessment	I am OK to commence works

Device Notes

Device Notes

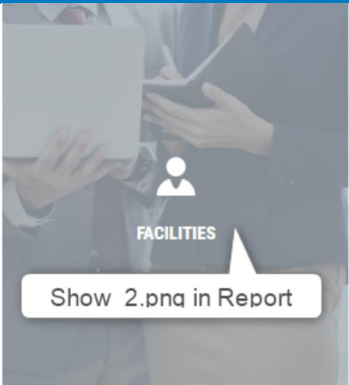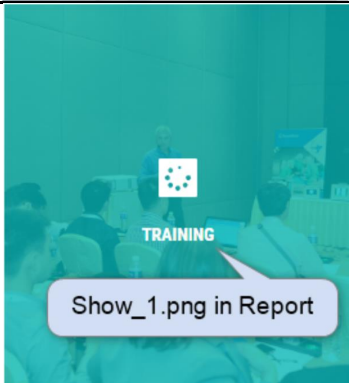

Image	Asset	Comments
	Air System: Oil Lub Screw: Medical: 7.5 Hp / 5.5 KW: Triplex: BeaconMedaes: 4233600592: mAIR-GTF4-722 HTM02-01 AC GA5 dMED35, s/n MNL_2498	Image Show_2.png uploaded
	Air System: Oil Lub Screw: Medical: 7.5 Hp / 5.5 KW: Triplex: BeaconMedaes: 4233600592: mAIR-GTF4-722 HTM02-01 AC GA5 dMED35, s/n MNL_2498	Image Show_1.png uploaded

Authorization to carry out work with/without an order number.

Notes: Customer will agree to supply Test Automatic UK with a valid order number or payment for non-account holders for the actual time spend, traveling and parts used during the breakdown intervention by next working day.

Customer Confirmation

Service Engineer Confirmation

Signed by Robot_UK
 (Being duly authorized to sign on behalf of the customer)

Signed by Robot_UK Tester



Test_Facility_Robot_UK
28-Jan-2025 005929

